

Kodak InSight


Dental Film in SURESOFT Packets

The first and only dental film with a welcome innovation: a specially designed cushioned edge for optimal comfort. It puts a soft, protective cushion between the edges of the film and your patients' gums. Now you can offer your patients the reduced radiation levels and enhanced safety of Kodak InSight dental film, with a new level of comfort. That's bound to make them smile.



The dental film with a cushioned edge.

The most comforting thing to happen in dentistry since the padded headrest.

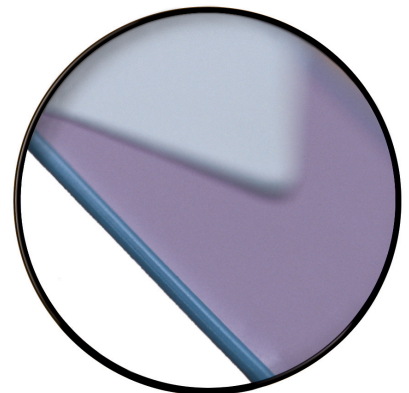

A new level of comfort for you and your patients

- A soft, round, cushioned edge covers the perimeter of the dental x-ray film packet.
- Reduces potential discomfort during x-ray exams and requires no special handling or preparation.

Plus, a combination of benefits that InSight film delivers

Short exposure time

InSight film is a true F-speed film for fast exposure times.

Up to 60% less radiation exposure

Count on enhanced safety for both patients and office personnel.

Superb image quality

Patented T-Grain® emulsion technology delivers images with exceptional sharpness and clarity.

Conversion is simple

From D-speed films, such as Kodak Ultra-speed film, reduce exposure time by 60%; from E-speed films by 20%. If you are currently using an amber safelight filter in your darkroom, change it to a Kodak GBX-2 safelight filter or upgrade to a new Kodak LED safelight.

Ordering Information

For KODAK ULTRA-SPEED Dental Film

Size	Description	Code	Films/Pkgs	CAT Number
Size 2	Kodak InSight dental film in SureSoft packets	IP-21S	80	811 8226

For more information, call **1-800-933-8031** or visit: www.kodakdental.com



Carestream Health, Inc.

© Carestream Health, Inc., 2008.
InSight, SureSoft, T-Grain and Ultra-speed are trademarks of Carestream Health, Inc.
The Kodak trademark and trade dress are used under license from Kodak.